

MR2=Manual Operation    MR3=Motorized or Manual



*MR2 & MR3 Supplied Complete with High/Low Profile Disposable Steel Blade Holder and Specimen Holder (choice of Cassette Clamp or Fixed Jaw Vise Clamp)*

*Wide range of accessories for Histology, Materials and Forensic Applications*


### FEATURES

- ◆ Heavy base casting to reduce vibration
- ◆ Vertical cutting stroke with solid cross roller bearing specimen carriage
- ◆ Digital Section Counter Display
- ◆ Touch Screen Control Panel with Easy to Clean Membrane Switches
- ◆ Section Thickness range 0.5 – 60 microns
- ◆ Electronic Feed with Fast Return
- ◆ Separate Trim Control Function with Footswitch
- ◆ Preset Trim Setting for Fast Block Trimming
- ◆ Removable Waste Tray


# SPECIFICATIONS

<b>Dimensions:</b>	H-279 mm x W-273 mm x L-457 mm Weight - 24.3 kg
<b>Maximum Specimen Size:</b>	55 x 50 mm
<b>Cutting Stroke:</b>	64 mm
<b>Retraction:</b>	70 microns
<b>Electronic Feed of Knife Stage:</b>	28 mm
<b>Handwheel Lock:</b>	Positive lock at 12:00 o'clock or 6:00 o'clock position by pulling out black hand wheel lock knob on left side of microtome.
<b>Specimen Orientation:</b>	x and y axis with 8 × tilt plus rotation
<b>Power Supply:</b>	Input 100-240 VAC, 50/60 Hz, 0.8 Amps Output 24 VDC, 36 Watts

## Manual

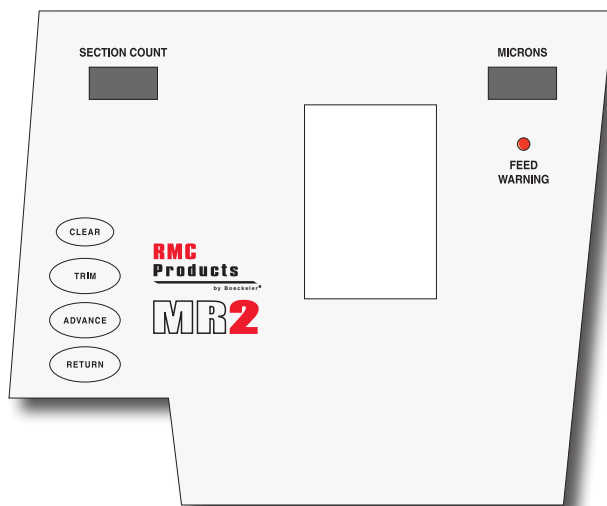

## Semi-Automatic

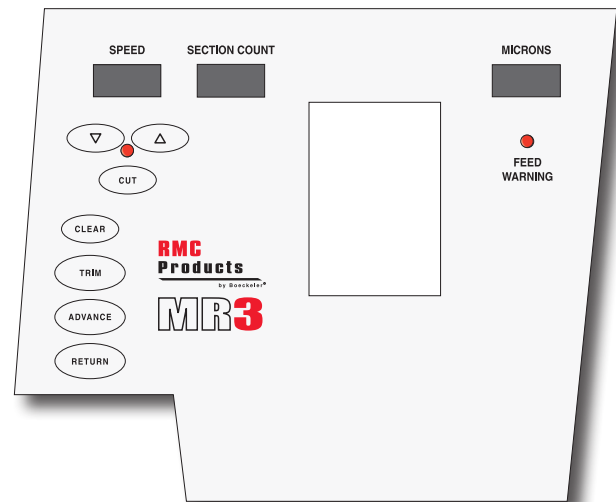

### Additional Specifications for MR3:

Variable cutting speed from 0 – 400 mm/second  
Footswitch control for both Cut & Trim Functions

Boeckeler Instruments, Inc.

4650 S. Butterfield Drive • Tucson, Arizona 85714 USA

Phone: 520.745.0001 • Fax: 520.745.0004 • E-mail: [info@boeckeler.com](mailto:info@boeckeler.com) • [www.boeckeler.com](http://www.boeckeler.com)

© 2008 Boeckeler Instruments, Inc. Boeckeler is a registered trademark of Boeckeler Instruments, Inc., Tucson, Arizona. All other trademarks are the property of their respective owners. All specifications subject to change without notice.